

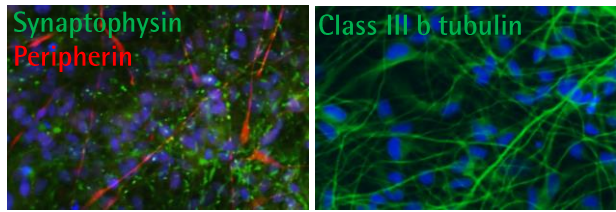


# A pre-clinical assay for testing the effects of compounds/drugs on innervated skin *ex vivo*

## 1. Differentiation of human iPSC derived neural stem cells *in vitro*


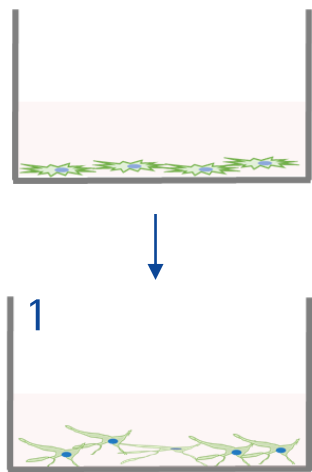


Human iPSC derived neural stem cells start to express **class III  $\beta$ -tubulin** which is associated with neuronal maturation (Sainath and Gallo, Cell Tissue Res 2015), **peripherin** which is a peripheral nervous system neuronal marker (Yuan et al., J Neurosci. 2012), and **synaptophysin** which is a marker of mature neurons (Kwon et al., Neuron. 2011)


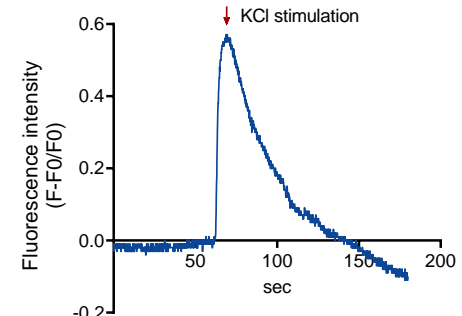


Differentiated human iPSC derived neural stem cells are functional. They release calcium in response to a depolarizing potassium chloride (KCl) stimulus. (response of n = 1 neuron is shown)



## 2. Initiation co-culture of human skin with human differentiated iPSC derived neural stem cells *ex vivo*

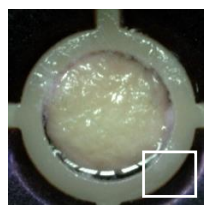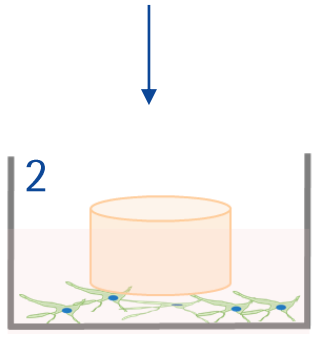



Topical application possible

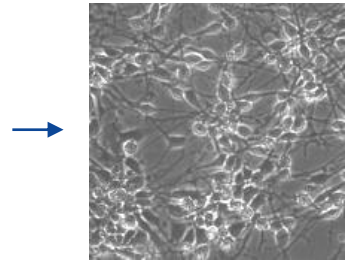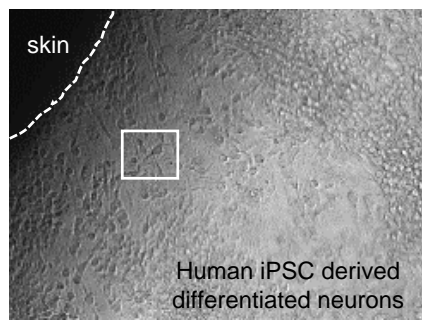

## 3. Model ready to use: Human skin punch is fully re-innervated

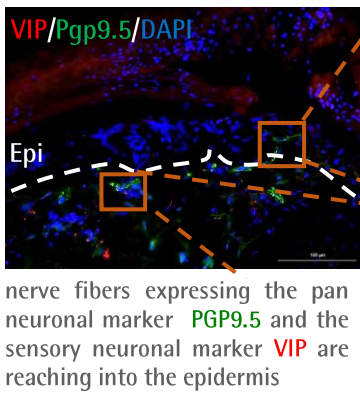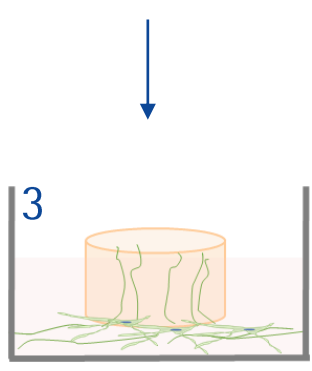

nerve fibers expressing the pan neuronal marker **PGP9.5** and the sensory neuronal marker **VIP** are reaching into the epidermis

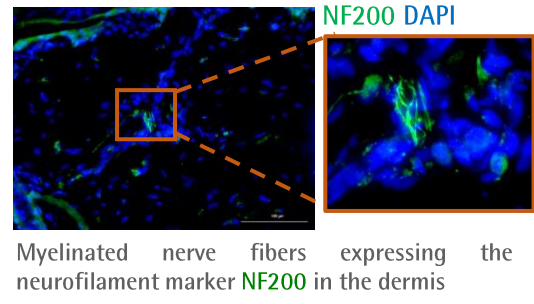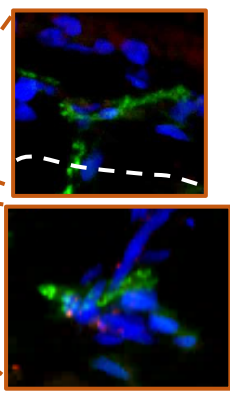

Myelinated nerve fibers expressing the neurofilament marker **NF200** in the dermis

This novel assay can be utilized for dissecting and manipulating the bi-directional communication between defined skin and hair follicle cell populations and (sensory) human nerve fibers under stringently controlled *ex vivo* conditions, or for testing cosmeceuticals or drugs that target the cross-talk between human skin and hair follicles and cutaneous nerve fibers.  
**Relevant for e.g. sensitive skin, itch, atopic dermatitis, psoriasis**

Keren et al., *Sci Adv.* 2022; Cheret et al., *JID* 2022

**Contact us for a customized study:**


Acting CEO:  
 Dr. Marta Bertolini (PhD)  
 CSO:  
 Dr. Janin Edelkamp (PhD)

m.bertolini@monasteriumlab.com  
 j.edelkamp@monasteriumlab.com  
 + 49 (0)251 93263-080

For more details see also our webpage:  
[www.monasteriumlab.com](http://www.monasteriumlab.com)

All about  
skin and hair  
bioscience!

State-of-the-art technology and expertise  
for all your pre-clinical, mechanistic, and  
clinical needs in dermatology research.

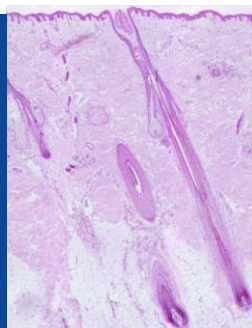

- ▶ Pre-clinical Research
- ▶ Clinical Research
- ▶ Education



"We combine  
our unique expertise,  
our project design creativity,  
and our passion to advance  
our clients' success in  
delivering novel and game-  
changing skin and hair  
research solutions"

Founder & CEO:  
Prof. Dr. Ralf Paus

Your one-stop  
source for all *in vitro*,  
*ex vivo* and  
*in vivo* testing  
and more.

Reasons why you should choose  
Monasterium Laboratory:

- Cutting edge methodologies and techniques
- Tailor-made & customized assays for all needs
- A focus on novel targets and therapeutics for skin & hair disorders: identify-characterize-validate
- Delivering mechanistic action stories, biomarkers & predictors of response
- Claims support for cosmetic ingredients in skin or hair follicle models
- Clinical trials carried out with strategic partners for healthy skin and hair benefits
- Comprehensive project reports & manuscript drafting and submission

**Monasterium Laboratory**  
Skin & Hair Research Solutions GmbH  
Mendelstr. 17, 48149 Münster, Germany

[www.monasteriumlab.com](http://www.monasteriumlab.com)

For enquiries, please contact:

Acting CEO:  
Dr. Marta Bertolini (PhD)  
[m.bertolini@monasteriumlab.com](mailto:m.bertolini@monasteriumlab.com)  
+ 49 (0)251 93263-080